

# HARVESTER 9600 Robot Accessible

From the Benchmark Family of 96-Well Harvesters

- » *Harvest to Filter Plates or Filter Mats (8 x 12 format on 9mm centers)*
- » *Fully Programmable*
- » *“Hot” and “Cold” Waste*
- » *Removable Wash Tips*
- » *Host or Stand-Alone Control*
- » *Proven Design*



## TOMTEC

## Filter Mats or Filter Plates

The Harvester 9600 is available in filter plate or filter mat formats for robotic operation. Standard glass fiber filter mats for the Betaplate® or MicroBeta® are carried in frames for robotic handling. UniFilter® plates for the Topcount® or MicroBeta® are handled directly. Special Hotels are available for either format.



Betaplate and MicroBeta are registered trademarks of Wallac Oy  
Unifilter and Topcount are registered trademarks of Packard Instrument Co.

## Proven Harvester Design

Tomtec's proven design for 96-well harvesters is incorporated into the Harvester 9600. It uses the same programmability and removable tip design for wash/aspirate and overflow control. Automation has simply been added to a very reliable design, that can be used in a stand-alone mode or controlled by an RS232 link with a host computer.

## Automated Buffer Supply and Waste



10-liter Buffer Bottle and Autotrap


Two 10-liter automated buffer supply bottles are offered. Each bottle has high and low level float control that interfaces with the Harvester 9600. These may be filled automatically from a remote source. Low air pressure (0 to 5 psi) is used to control buffer flow rate to the Harvester 9600.

Two solenoid operated Autotraps complete the waste collection. The Autotraps are float controlled and can switch from the "trap" to "drain" function as the program requires.

## Ordering Information

- Harvester 9600 MachIII  
8 x 12 MicroBeta® Format  
(1/60/115vac)  
Cat. No. 696-300
- Harvester 9600 MachIII  
UniFilter® Plate Format  
(1/60/115vac)  
Cat. No. 696-300U
- Remote Operator Panel  
Cat. No. 696-325
- Reagent Container with fittings  
(10-liter)  
Cat. No. 696-550-1
- DC Autotrap with ½" ID tubing,  
24 volts (from Harvester 9600)  
Cat. No. 96-15-1

## Dimensions in Inches

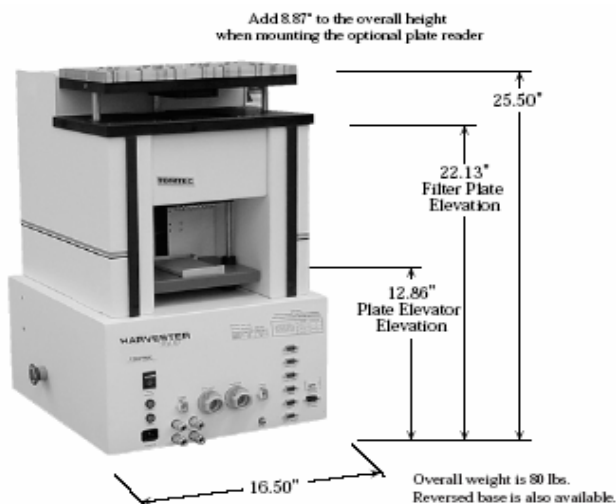

## TOMTEC

1000 Sherman Avenue, Hamden, CT 06514  
203.281-6790 Fax 203.248-5724  
Toll Free: 877-TOMTEC3 www.tomtec.com